



UnicLine
Mobile
- Real
Mobility

UnicLine Mobile
- you will not
feel the
difference to
stationary
equipment



Ideal for dental care at
schools, nursing homes,
companies and army
facilities as well as for
emergency treatment.
Practically all places
where there is a need for
mobile treatment without
having to compromise.



UnicLine mobile
- from boxes to
clinic in less than
15 minutes

Same hygienic and
ergonomic conditions
as in a stationary
clinic.



UnicLine[®]
- FEEL THE DIFFERENCE



Multi armrest

The UnicLine patient chair can be delivered with one or two armrests.



Joystick

Individual placement of up to 3 joy-sticks offers you and your assistant to maintain ideal working postures.



Easy to clean

The UnicLine upholstery is easy to dismount for cleaning.



Ergonomics

Arm in cast aluminium with metal springs for easy handling without use of brake button or the like.

Instrument table with 2 micromotors, 3-function syringe, ultrasonic scaler and LED light polymerization lamp.


Display with readings of speed, direction of rotation, choice of spray, clock and timer operation.



Trays can be mounted on both sides.

3 selective suction hoses suspended in a flexible arm.

Perfectly balanced suspension of instruments, with extra long hoses.

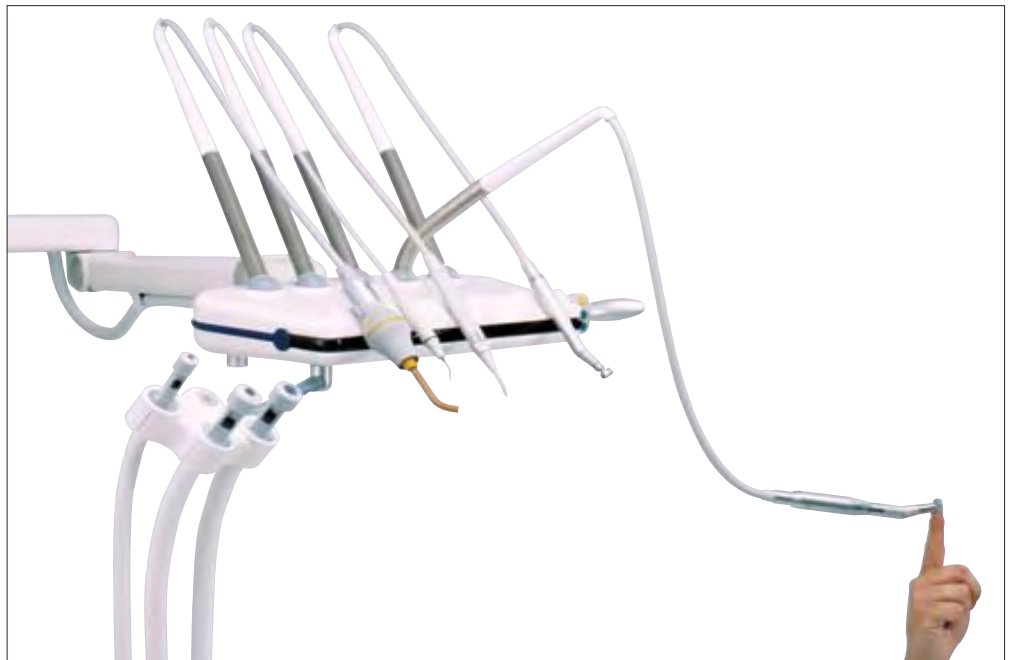



Everything on wheels

The transport boxes are specially designed with separate spaces lined with polyester foam.

All parts are stored firmly and safely during transport. Arm, instrument table, foot control and OP Lamp in one box. X-ray with timer in separate transport box.



UnicLine[®]
- FEEL THE DIFFERENCE

Silenced unit cabinet containing suction motor with amalgam separator, compressor, outlet container with odour filters and a water container.

Specifications UnicLine Mobile

Micromotors – 2 pcs Bien Air MC3 with light
Rpm 40,000
Torque 2.9 Ncm
Bulb 3.5 VDC; 0.75 A
Length 94.3 mm including attachment
Weight 98 gm

3-function syringe – Luzzani
Weight 80 gm
Autoclavable tip and casing

Ultrasound scaler – EMS Piezon
Including 3 tips
Length 104 mm
Weight 40 gm

Curing light -LED-polymerization lamp mectron starlight S
Weight only 75 g
5W LED with patented optical focusing
Light intensity Ø 8 mm light guide > 1.000 mW/cm²
Wavelength range 440-480 nm with peak at 460 nm
Working life of LED 2.500.000 cycles, 10 seconds each

OP Lamp – Faro Forma 2000
27,000 lux measured at 700 mm from surface
Colour temperature 5000°K
Bulb power 95 W

Suction motor –Dürr VSA300
With connections for 3 selective suction hoses,
2 pcs of 15 mm, 1 pc of 11 mm
Integrated amalgam separator
Capacity 300 l/min
Max liquid flow 9 l/min

Patient chair
2 spindle motors
Max weight on the chair 240 kg
Max lifting height 920 mm
Min height 510 mm
Weight of chair 95 kg
Including transport wheels and lifting bar

X-ray – Satelec X-mind A/C (From medio 2004)
Classification Electromedical equipment, class 1 type B

High voltage 70 kV
Supply voltage 220/230/240 V – monophas 50/60 Hz
Focal spot 0.7 mm

Transport box – Arm system
2 transport wheels
4 small wheels for short removals
8 lifting devices
2 supporting legs for upright storage
Protection edgings of aluminium
Dimensions L x B x H – 108 x 71 x 40 cm
Weight with contents 80 kg

Transport box – X-ray
2 transport wheels
5 lifting devices
Box and lid are lined with polyester foam
Protection edgings of aluminium
Dimensions L x B x H – 106 x 33 x 27 cm
Weight with contents 48 kg

Compressor – Jun-Air with microfilter
Capacity 50 l/min
Max pressure 8 bar
Tank 4 l
Noise Level 45 db

Outlet container – 2.5 l container
With sound warning when the container is almost full and
automatic cut-out when the container is completely full.

Water container – 1 l container
with quick coupling

Unit cabinet
2 supporting legs
2 transport wheels
Transport handle
Dimensions L x B x H – 60 x 58 x 87 cm
Weight with contents 86 kg

Total dimensions
Weight 309 kg
Volume (packed) 1.4 m³
Space needed (mounted) 5 m²

heka*dental*[®]



Heka Dental A/S
Baldershøj 38
DK-2635 Ishøj, Denmark
Phone +45 4332 0990
Fax +45 4332 0980
mail@heka-dental.dk
www.heka-dental.com

Future

For further information - please contact your depot,
service technicians or Heka Dental.

With reservation of
technical changes and
colour changes.